


# Reusable EEG Cup Electrode Silver/Silver Chloride

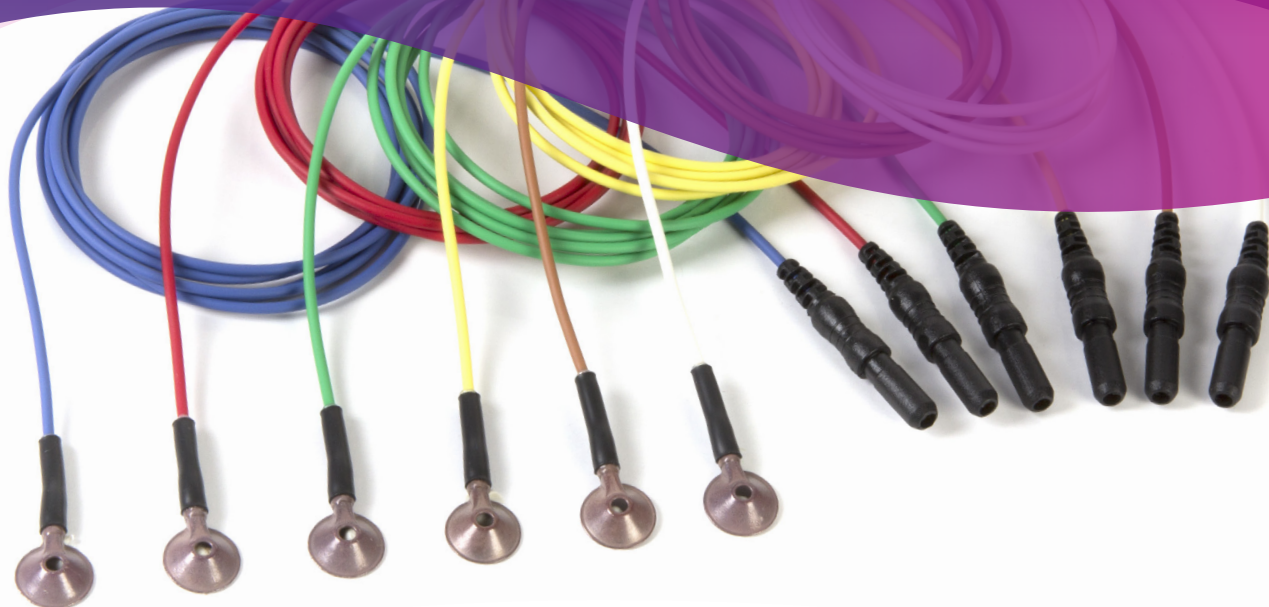

**At Technomed we understand the importance of a crystal clear signal, ease of use and electrode durability.**

**Our new Reusable Ag/AgCl EEG Cup Electrode is designed to give you reliable signal, and ease of use time after time.**

**The new electrode head is designed with a larger dome volume for longer term use, and the connector is optimized for easy headbox insertion and extraction.**

#### **User Benefits**

- Clear and reliable signal with optimized layer of frequency stable Ag/AgCl on a solid Silver cup
- Optimal patient comfort with soft and flexible lead wires
- Extra large dome volume for easy filling of cups and longer term use
- Optimized connector for easy headbox insertion and extraction
- Long electrode life with Kevlar reinforced wires

## Available variants

Materials	Wire Length	Part Number
Silver/Silver Chloride (Ag/AgCl)	1 m	TE/C22-434/99
Silver/Silver Chloride (Ag/AgCl)	1.5 m	TE/C22-634/99
Silver/Silver Chloride (Ag/AgCl)	2 m	TE/C22-834/99
Silver/Silver Chloride (Ag/AgCl)	2.5 m	TE/C22-934/99

All variants are packaged in sets of 12, 2 x 6 colors: red, green, white, yellow, blue, brown



## Specifications

### Dimensions

Diameter of cup	10 mm
Height of cup	3.3 mm

### Properties

Lead-free	√
Latex-free	√
BPA-free lead wire	√
DEHP-free lead wire	√
RoHS compliant	√
Lead wire general	Tinsel wire
Lead wire core	Kevlar and Silver plated Copper
Lead wire insulation	TPU
Spring loaded connector	Gold plated Brass/ PVC
Connector	DIN 42802

### Cleaning and disinfection

Cleaning by washing with non-scratching tools and lukewarm water. Disinfecting with alcohol.

### Precautions



Your Technomed distributor: